

Low Maintenance

The permanent magnet does not require cryogenics, weekly servicing, or preventative maintenance. With standard 5mm NMR tubes, and automated shimming routines, the 100 can be plugged into standard wall power and used daily without requiring trained staff to maintain it.



Superior Performance

The Nanalysis 100 MHz offers 25% higher resolution than competing products. This affords better peak separation and superior sensitivity, while still maintaining an accessible, all-in-one, easy-to-use footprint.



Compact with Superior Ergonomics

The 100 MHz is the next evolution in benchtop NMR spectrometers. It can be located on any bench with access to a wall plug-in. The tilt and swivel display provides optimal touch screen access from any angle, whether sitting or standing. Easy access connection ports in both the front and rear of the spectrometer allow for easy data transfer access.



Easy Integration

The 100 MHz has been designed specifically to be incorporated into your existing laboratory. With numerous secure and accessible connectivity options, data can be easily exported and worked up in third party software programs.

- Implementation of optional software IQ/OQ
- Optional software packages allow for restricted user access
- Incorporation into LIMS
- Individual user logins
- Project creation and even simplistic design and implementation of applications

Inquire for more information about these software packages.



Unlike any benchtop NMR you've seen before...

With the same ease of use, and all-in-one convenience that Nanalysis has pioneered, the 100e brings a new level of power right to your benchtop. The 100e, everything you need to bring powerful NMR right to your lab.

Technical Specifications

Operating Frequency:	100 MHz (2.35 T)
Magnet:	Permanent, no cryogenics
User Interface:	Integrated computer and/or remote accessibility
Nuclei:	^1H or ^{19}F
Lock:	Internal ^1H or ^2H
Sample:	Standard 5mm NMR tubes
Compatibility:	JCAMP-DX, Mnova, Delta, ACD/Labs, Topspin, Labview, Matlab
Resolution & Lineshape:	<1.0 Hz
Sensitivity:	>220 : 1 (1% EtBz)
Operating Temperature:	16 to 28°C
Power Supply:	100-240 VAC, 50-60 Hz

Dimensions:	14.6 x 16.3 x 25.8" (plus screen) (37.1 x 41.4 x 65.4 cm)
Weight:	215 lbs (84 kg)

Nanalysis Corp.
Bay 1, 4600 5 Street NE
Calgary, Alberta, Canada T2E 7C3

+1.403.769.9499
Sales: sales@nanalysis.com
www.nanalysis.com



 **nanalysis**
100e




Customizable Experience

Pre-loaded with many standard, fully-adjustable experiments, users can immediately acquire results, and customize parameters. Monitor and assess reaction completeness, purity and relative composition right in the lab.


Versatile User Interface

The instrument can be controlled via the touchscreen, with a keyboard and mouse, or through an external computer to:


- **Load Standard Experiments**
The standard experiments are (1D, COSY, JRES, T_1 , T_2 , etc.)
- **Design Pulse Experiments**
Our advanced software allows users to easily create and modify pulse programs, including pulse power, timing, phase cycling, etc..
- **Queue Experiments**
Through automated analysis, set up a series of experiments to be acquired automatically (e.g. 1D, COSY, JRES, T_1)

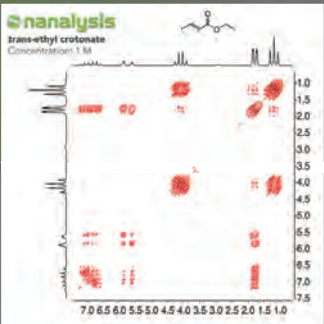


COSY

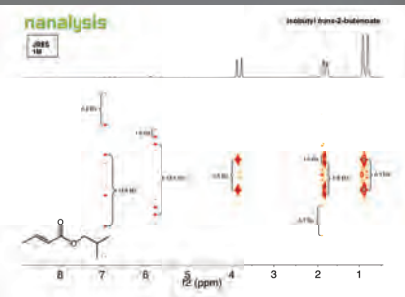

JRES

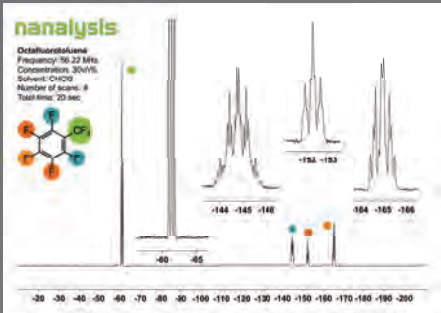

19F

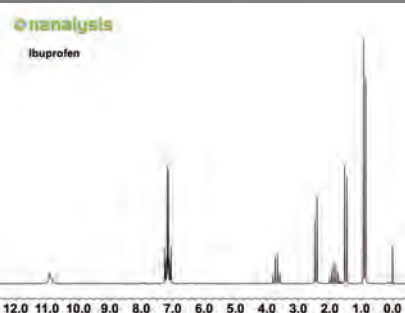

1H

Quick and Easy-to-use

Like the rest of the NMReady family, the integrated touchscreen permits users to load standard experiments directly or modify key acquisition parameters. Data can be worked up on board with the included processing software, easily exported for use with third party programs, or printed. With superior sensitivity, samples can be collected in seconds.

1 Configure

Experiments are selected from simple dropdown menus and loaded with standard acquisition parameters. These can be modified as desired from a series of convenient pop-up menus.

2 Acquire

Once configured, data is collected by merely selecting 'Go'. Each scan is collected in seconds, depending on acquisition parameters. The first scan is displayed and updated every additional four scans.

3 Analyze


State of the art on-board data processing makes it easy to process and analyze 1D, 2D spectra on the instrument. Additionally, spectra are saved as standard JCAMP-DX files and are compatible with many third-party software tools.

Connectivity

Ethernet / WiFi, USB, Serial, HDMI



HDMI

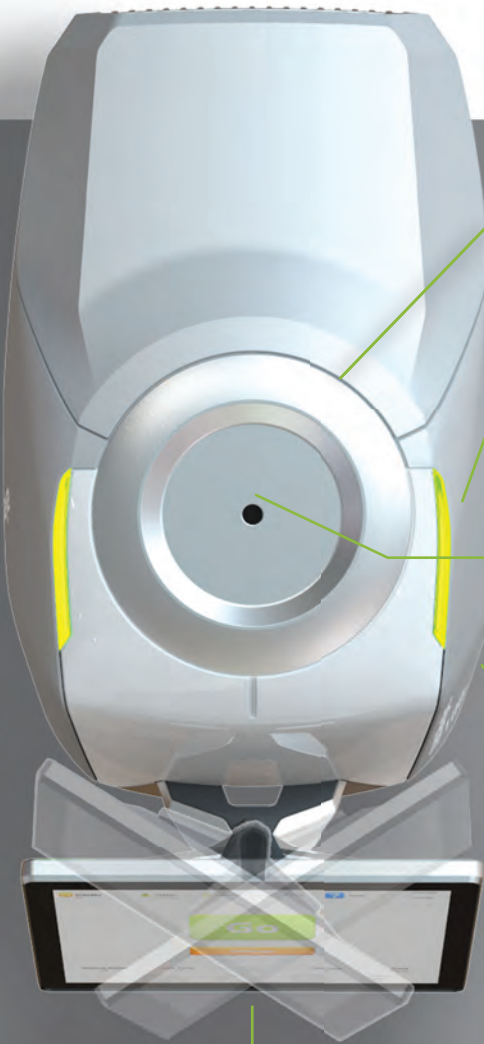

Innovative magnet design
<2G outside

Progress Indicator
Large status light so you know when a scan is finished from across the room

Sample Access Port
5mm NMR tubes

Enhanced stability
Vibration absorbing feet and precision machined structures for optimal stability

Ergonomic Display
State of the art multiposition possibilities

Quick Access
Front USB ports and power button



Optimal Thermal Control
External heat sink and side vents for built in convection and conduction heat management

